

Operation Instructions

Utility Tray (9A96001) used on Podiatry Tables 414, 416, and 417

1. Position the podiatry table to the full tilt position.
2. Disconnect the podiatry table's electrical supply cord to prevent electrical shock.



WARNING




EQUIPMENT ALERT

Failure to disconnect the podiatry table's electrical supply cord could result in personal injury or equipment damage.

3. Find the two holes used to fasten the bracket weldment (1) to the table.
4. Position the bracket weldment.
5. Securely fasten the two 1/4" bolts (2) supplied.
6. Adjust the tray assembly (3) to the desired position.

NOTE

Tray assembly can be rotated 90 degrees.

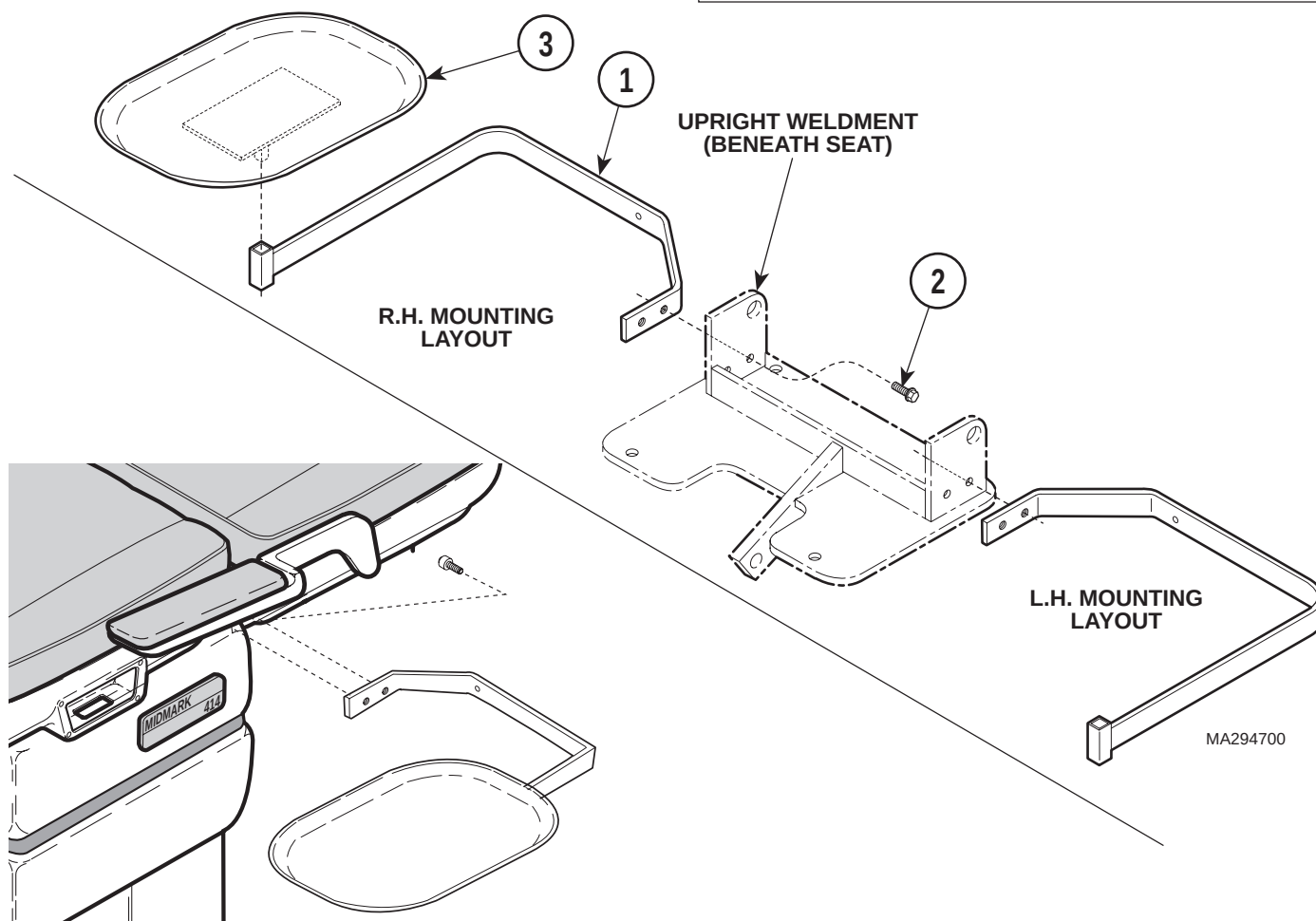

MA31060

MIDMARK®

Midmark Corporation, Versailles, Ohio 45380 U.S.A.
937-526-3662 TELEX 288007 FAX 937-526-5542